



WHILL Model C Lite

The World's Only Lightweight, Foldable Carbon Fiber Wheelchair with Omni-Wheels



Lightweight.
High Performance.
Comfortable.
Adjustable.



Cutting Edge Features



Incredibly Lightweight

The world's only foldable electric wheelchair with a molded carbon fiber frame. It is the lightest in its class.




Best Terrain Coverage

Powerful motors combined with large front wheels deliver the best terrain coverage in the lightweight category.



Fully Adjustable

Adjustable back support angle and height, armrest height, plus foot support height and length.

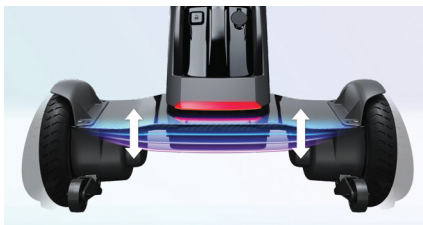

Carbon Blade Suspension

Flexible blade absorb shocks and vibrations effectively, ensuring a smooth ride even on solid tires.

The World's Only
Lightweight, Foldable
Carbon Fiber Wheelchair
with Omni-Wheels





Quick Folding: Fold in just 2 seconds. Once folded, it stands upright for easy storage and transport and you can roll it effortlessly.



Easy One-Hand Folding: A new folding mechanism lets you fold with one hand. No awkward two-handed bending required.



Stay Straight on Slopes: Sensors independently adjust each rear wheel, keeping you moving in a straight line on slopes without drifting.




Maneuverability: A 31.9-inch turning radius makes tight turns in narrow spaces effortless.



High Obstacle Clearance: Clear curbs up to 2-inch with ease thanks to large front wheels and two powerful motors.



Largest Legroom in Its Class: Flip-up footrests and generous foot support improve usability for a wider range of users.

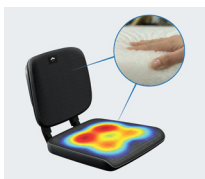

Medical-Grade Cushion: The same premium material found in Japan's bullet train first-class seats delivers all-day comfort.




IPX5 Water Resistant: Rides confidently in rain and outdoor conditions with IPX5-rated water resistance.

Get the Most Out of Your WHILL



Load Into Your Car with Ease

The optional Loading Arm give you extra leverage when lifting the device into your vehicle.

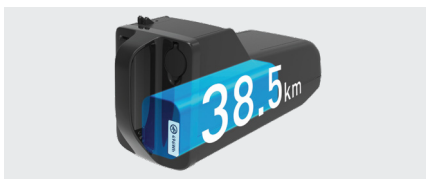

Long-Range USB Type-C Battery

⚡ 38.5 km (23.9 miles) driving range

⚡ Developed in partnership with BYD

100% Recycled Carbon Fiber

Made with recycled carbon resin, reducing CO2 emissions.



Customization



○ Iconic White



● Black



● Light Green



Basket



**Long-Range USB
Type-C battery**



High Back Support



Loading Arm

Specifications

Driving Range with Full Battery Charge *1	Up to 23.4 km (14.5 miles)
— with optional Long Range Battery *1	Up to 38.5 km (23.9 miles)
Maximum Speed	6 km/h (3.7 mph)
Ground Clearance	3" (75 mm)
Obstacle Clearance	2" (50 mm)
Maximum Incline	10 degrees
Turning Radius	31.9" (811 mm)
Weight Capacity	300 pounds (136 kg)
Air Travel Approved	Yes
Weather Tested	IPX5
Suspension	4 Wheel Independent
Seat Size	W 16.5 x D 18" (W 420 x D 450 mm)
Seat Height	20" (508 mm)
Device Width	23" (587 mm)
Device Length	38.5" (978 mm)

Device Height	31.3" to 34.2" (795 to 868 mm)
Device Weight *2	41 lbs (18.6 kg)
Folded Device Width	23" (587 mm)
Folded Device Length	20" (500 mm)
Folded Device Height	30" (755 mm)
Battery Weight	6.6 pounds (3 kg)
Charging Time	Up to 4.5 hours
Battery Type	Lithium-ion (25.6 V / 10.6 Ah)

*1 The driving range of this device varies depending on the types of roads or terrain it is driven on, as well as the load, ambient temperature, operating method, and other factors.

*2 Without the battery

Warranty:
Frame: 5 years/ Seat and back cushion: 2 years/
Electrical and mechanical components, batteries,
and tires: 1 year. For complete warranty terms
and conditions, visit whill.inc/ca/warranty-policy

Register your device online at www.whill.inc or scan QR code.

The type class of the powered wheelchair is Class 1.



WHILL Inc.

410-433 Main Street, Winnipeg, MB, Canada R3B 1B3
(844) 699-4455 ca.info@whill.inc www.whill.inc



740-00247

The intended purpose of the WHILL Model C Lite powered wheelchair is to provide outdoor and indoor mobility to persons limited to a seated position that are capable of operating a powered wheelchair.

Do not use the device as a seat in an automobile or other vehicle. There is a risk of damage to the device, causing an accident.